# **MicroCLAVE®**

# NEUTRAL DISPLACEMENT CONNECTOR A Needlefree Closed System Device

Despite its diminutive size the MicroCLAVE Neutral Displacement Connector can be one of your most powerful allies in the fight against catheter-related bloodstream infections. Neutral displacement and saline flush help promote better clinical outcomes.



# Small is Powerful!

**SALINE***FLUSH* 

## **MICROBIAL BARRIER**

A dedicated internal fluid path makes the CLAVE Technology the only needlefree system proven to reduce catheter colonization in a randomized clinical trial (Bouza, 2003). CLAVE has a proven track record for assisting in the fight against CRBSI (Maragakis, 2006).

# **CATHETER PATENCY**


Neutral displacement virtually eliminates blood reflux into the catheter lumen upon disconnection of an administration device. No blood reflux is what makes the use of a saline flush possible.

## **PROTOCOL NEUTRAL**

The MicroCLAVE Neutral Displacement Connector supports many styles of clinical practice, eliminating the educational burden and risk of error due to variations in technique.

# Simple. Safe. Secure.

# MicroCLAVE® NEUTRAL DISPLACEMENT CONNECTOR

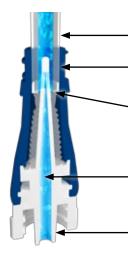

#### Administration Syringe

Outside luer of syringe is never in contact with the fluid path

Dynamic seal between tip of syringe and split-septum close the fluid path

Dedicated internal fluid path of the blunt cannula

Out to catheter hub

# **Small is Powerful!**



The small size, limited deadspace, and high flow rate of the MicroCLAVE Neutral Displacement Connector make it ideal for use in all areas of vascular access:

- CVC
- PICC
- Implanted Ports
- PIC
- Arterial



#### **TECHNICAL SPECIFICATIONS**

JE

Residual/Deadspace Volume: 0.04mL Flow Rate at Gravity: 165mL/minute Functional Activations: 600 Blood Compatibility: Yes MRI Compatibility: Yes Drug Compatibility: Alcohol, Lipids, Chlorhexidine, Chemotherapy

#### **ORDERING INFORMATION**

The MicroCLAVE is available in a wide variety of configurations to meet your facility's needs and requirements.

For more information on why the MicroCLAVE is right for you, contact your ICU Medical representative at 800-824-7890 or online at www.icumed.com.

U.S. Put. Nos; 5,685,866; 5,690,612; 5,694,686; 5,695,466; 5,873,862; 5,901,942; 5,928,204; 6,132,403; 6,572,592; 6,599,273; 5,738,663; 6,019,748; 6,325,782; 5,700,248; 6,132,404; 5,199,947; 5,344,414; 5,688,254; 5,776,116; 5,330,450; 5,971,950; 6,290,688; 6,669,673; 6,635,044; 6,572,592; 6,682,509; 6,726,672; 6,758,833 Other Patients Pending M1-1113 Rev. 05